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I T A L I A

Premium Features	3
Safety Precautions	3-4
Chair Components	4
Massaging positions	4
Chair Installation & Moving	5-6
Grounding Instructions	7
Chair Specifications	7
Operation Guide	7
Massage chair control buttons	8-11
Troubleshooting & Product Maintenance	11
Warnings	12
Other ANATOMICO products	13-14

Premium features

Therapeutic massage includes: Shiatsu (acupressure), Swedish, kneading, rolling, tapping, clapping, vibration, dual-action and tri-action techniques.

- Innovative auto scanning technology in backrest automatically or manually customizes a personal massage specifically for each user.
- Invigorating air pressure system in back, seat, calves and feet provides a total body massage.
- Six pre-programmed massage sessions with customizable strength, speed and intensity.
- Ability to customize massage by controlling speed, adjusting rollers and choosing partial massage areas.
- Three partial massage options in upper, middle and lower back.
- Shiatsu (acupressure) massage functions in seat, calves and feet.
- Narrow or wide back massage adjustment.
- Extended footrest adjusts to fit user's height.
- Automatic angle adjustment in backrest and footrest.
- Programmable timer with 15-30 minute automatic shut-off and pause mode.
- User-friendly LCD remote controller with graphic interface.
- Ergonomically designed for optimum comfort.


Safety precautions

When using the massage chair, basic precautions should always be followed to reduce the risk of electric shock, burns, fire or injury to persons. Please read the following precautions before using the massage chair.


- Make sure the power plug is securely inserted into a proper electrical outlet to reduce the risk of short-circuit and fire.
- Always unplug the massage chair from the electrical outlet immediately after using and before cleaning.
- DO NOT use the massage chair in areas where the humidity is high such as a sauna or steam room.
- DO NOT use the massage chair outdoors and avoid areas with direct sunlight.
- If the electricity is off, turn the power switch to the OFF position and unplug the massage chair from the electrical outlet to prevent any damages to your chair when the electricity comes back on.
- Never plug/unplug the power cord plug socket with wet hands.
- Close and continuous supervision is necessary when the massage chair is used by, on or near children or disabled persons.
- No more than one person should use the massage chair at a time.
- The massage chair is intended for use as described in the Instruction Manual. Do not use any attachments other than those recommended by the manufacturer.
- DO NOT attempt to stand or sit on the backrest, armrest or footrest in order to prevent accident, injury or damage to the chair.
- Avoid falling asleep while using the massage chair.
- Always keep the massage chair clean, and never insert any objects into or around the gaps of the massage chair.
- DO NOT move the massage chair by pulling on the power supply cord.

- Make sure there are no children, pets or other obstacles behind or underneath the massage chair.
- If you find any rips or tears in the fabric of the massage chair contact the service department and stop using the chair immediately.
- This massage chair is designed for household use and should not be used in non-residential settings.
- Please consult your doctor if you are in doubt as to whether or not you should use this massage chair.
- Please consult a doctor before using the massage chair if you are pregnant, undergoing any medical treatment, using a medical electronic device such as a pacemaker, or are experiencing any of the following: Malignant tumor, heart disease, back pain, abnormal or curved back bone, osteoporosis, or acute illness.
- If you experience pain while using the massage chair, stop using the chair immediately and consult your doctor.
- Never use a massage function on swollen or inflamed parts of the body.
- DO NOT use the massage chair for more than 30 minutes at a time.
- DO NOT use the massage chair where aerosol (spray) products are being used or where oxygen is being administered.
- This massage chair is a non-professional product designed to provide a comfortable massage and should not be substituted for appropriate medical treatment.

Chair components

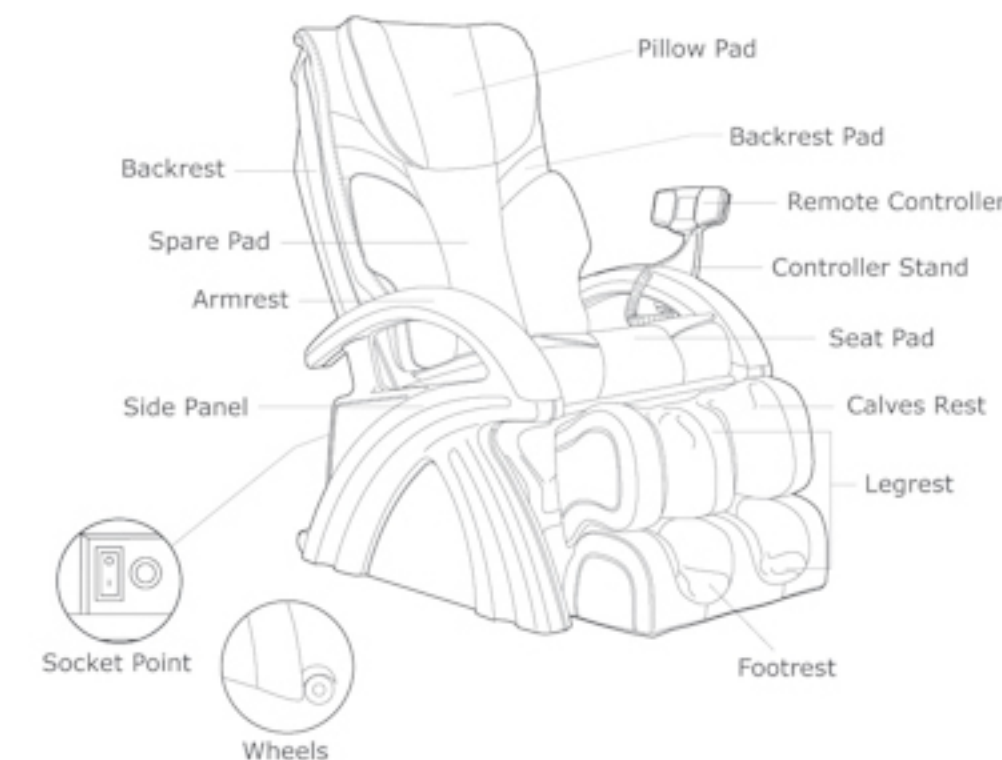
Massaging position

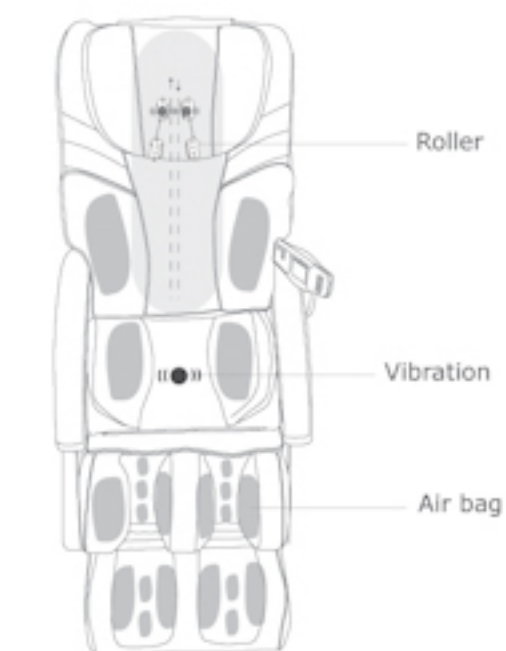
Chair installation

Step 1: When positioning the chair, make sure there is at least three feet of space behind the backrest and three feet of space in front of the footrest.

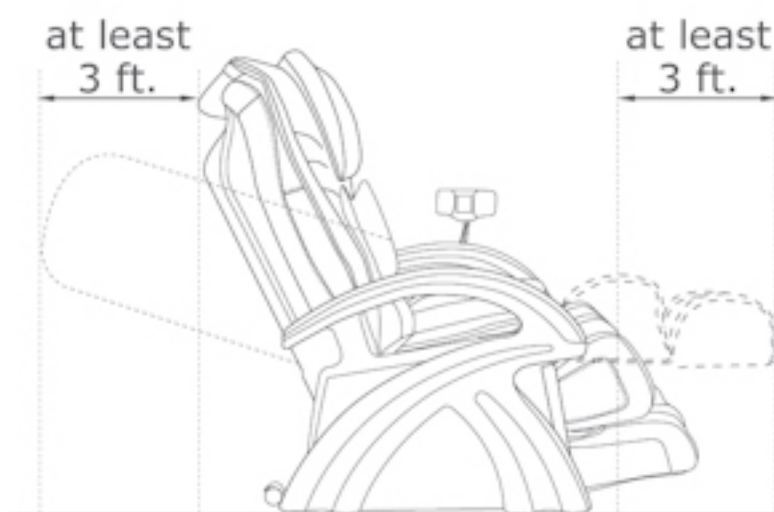Step 2: Lift the backrest, then lower the backrest into the seat brackets, Securely lock the flat head bolts into sides of the seat.


Note: Carefully unwrap the armrest to avoid scratching during installation.

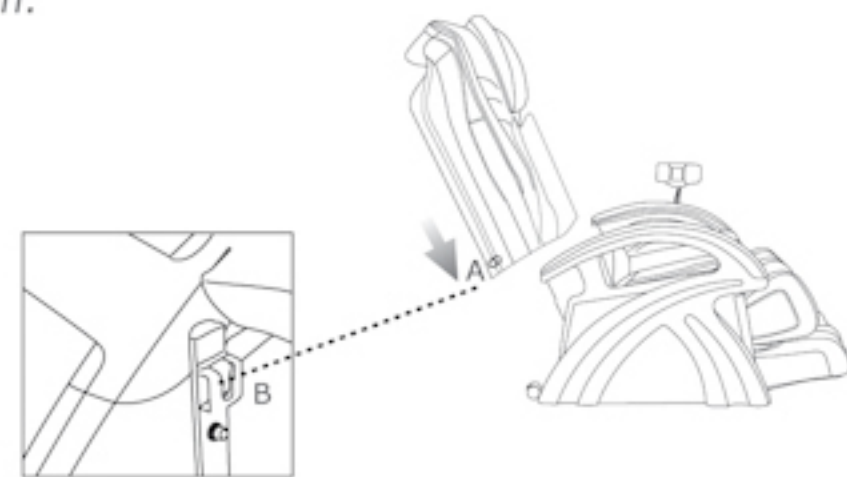A: Flat Head Bolt (on the backrest)
B: Holder (on the seat part)


Attach the connecting rod and the U stay fork of the reclining actuator with the flat head rivet, and secure with a split pin as demonstrated in PIC 1 and PIC 2.

- A. Connecting rod (on the backrest)
- B. U stay fork of reclining actuator (on the seat)
- C. Flat head rivet
- D. Split pin

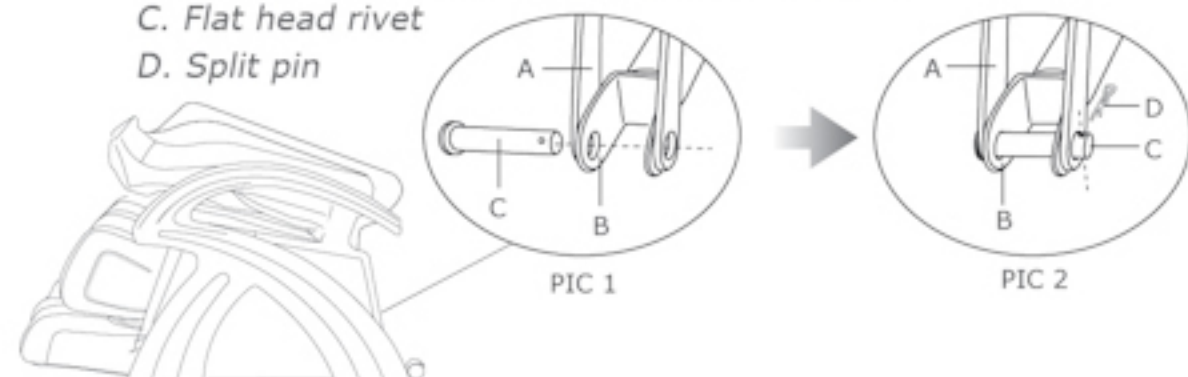
Step 3: Insert the wire plugs into the corresponding place in the seat according to the marks given.

⚠ Insert black to black and gray to gray, otherwise damage to the chair may result.

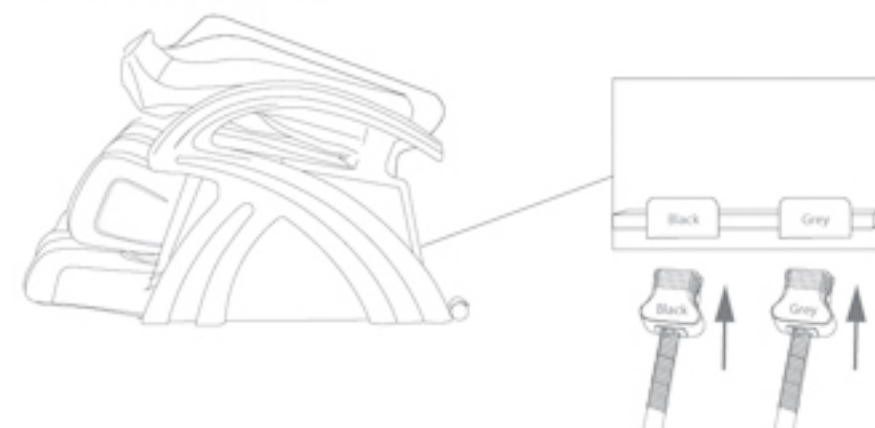Step 4: Attach the backrest pad and the pillow pad with the zipper.

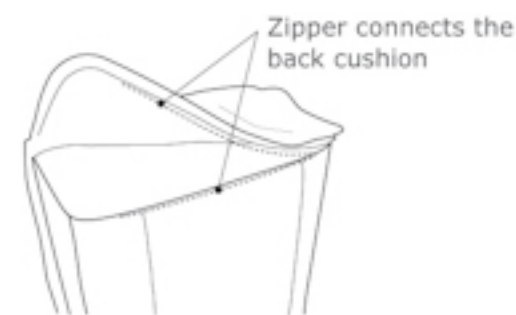
Step 5: Install the controller stand.

Tighten the controller stand screws.

Step 6: (optional) Remove the back massage strength-control pad from the backrest cover for a deeper massage.

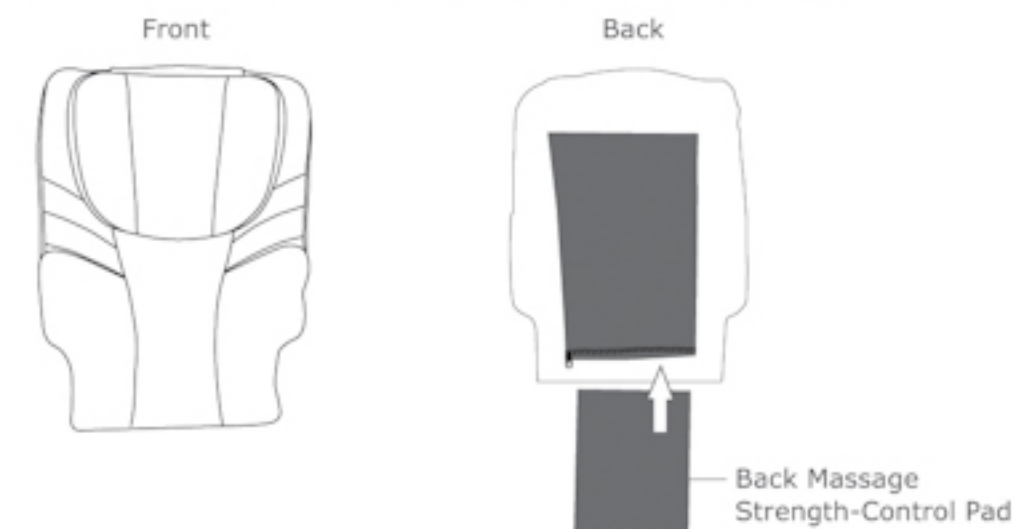
Note: The spare pad has been inserted into the backrest cover. If strong massage is desired, remove the spare pad from the cover, this strengthens the massage. If milder massages are desired, insert the spare pad into the cover; this lessens the strength of the massage.

Moving the chair

1. Restore the footrest to its original position.
2. Recline the chair to tilt approximately 45 degrees. Move the chair using the wheels as the fulcrum.

Note: 1. Always return the footrest to its original length and restore it to its original position before moving the chair. Otherwise, injury during movement may result.

2. Do not move the chair while a person is seated on it.
3. Place a mat on the floor and gently roll the chair over it to avoid damaging the floor.

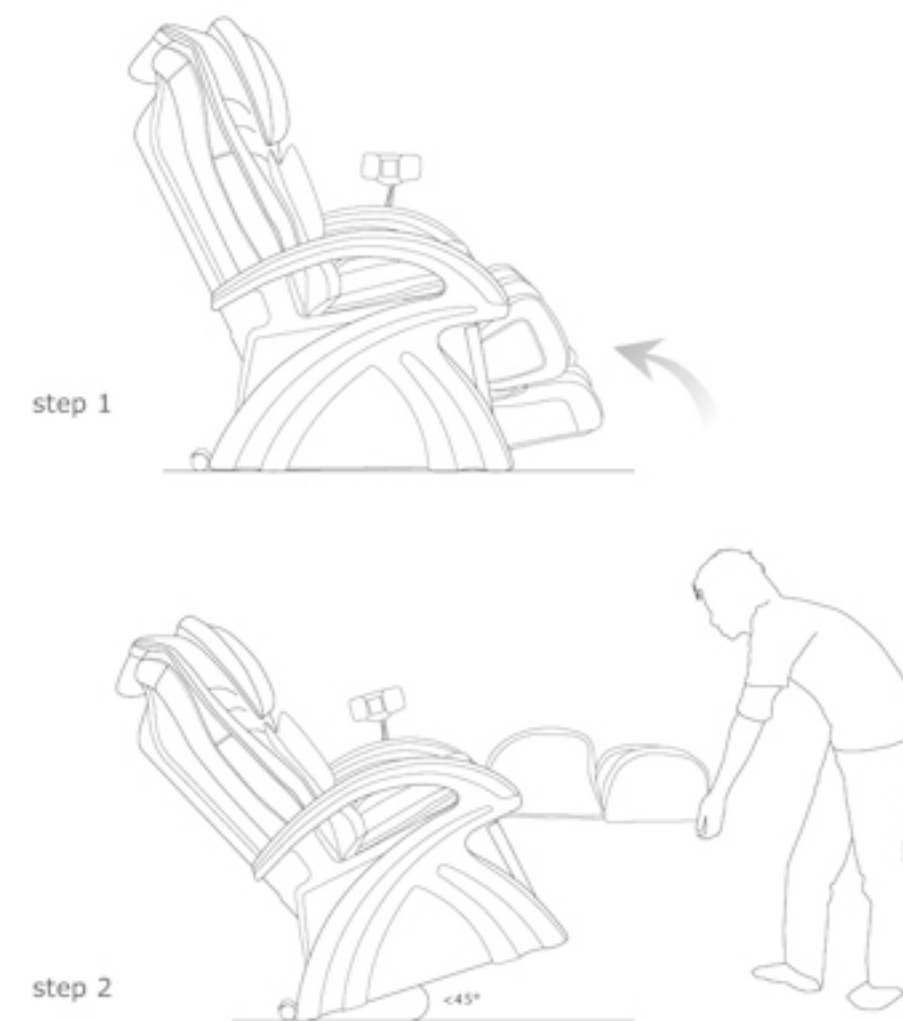
Grounding instructions

Fully insert the plug into an appropriate outlet. An improper equipment-grounding conductor may result in electric shock. Contact a qualified electrician if you are not sure that the product is properly grounded. Do not modify the plug that is included with the massage chair. If the plug does not fit the outlet, contact a qualified electrician to install an outlet to fit the massage chair plug.

Chair specifications

Description	Specifications
Rated Voltage	120V AC
Rated Frequency	60Hz
Rater Power Consumption	250W
Timer	5, 10, 15, 20, 25, 30 minutes
Max Time Rate	60 minutes
Dimensions(LxWxH)	Shipping dimensions (seat): 48" x 30" x 29"
	Shipping dimensions(Backrest): 44.5" x 25.5" x 16"
	Upright dimensions: 79.5" x 28" x 29.5"
	Reclined dimensions: 56" x 28" x 46"
Weight	Net weight (seat): 132 lbs
	Net weight (backrest): 51 lbs
	Gross weight (seat): 186lbs
	Gross weight (backrest): 62 lbs
Length of Wire	Controller wire: 1.3 m
	Power supply wire: 2.2 m
Usage Condition	Temperature: 10°C~40°C
	Humidity: 30~85 RH
Storage Condition	Temperature: 20°C~60°C
	Humidity: 30~85 RH

Operation guide

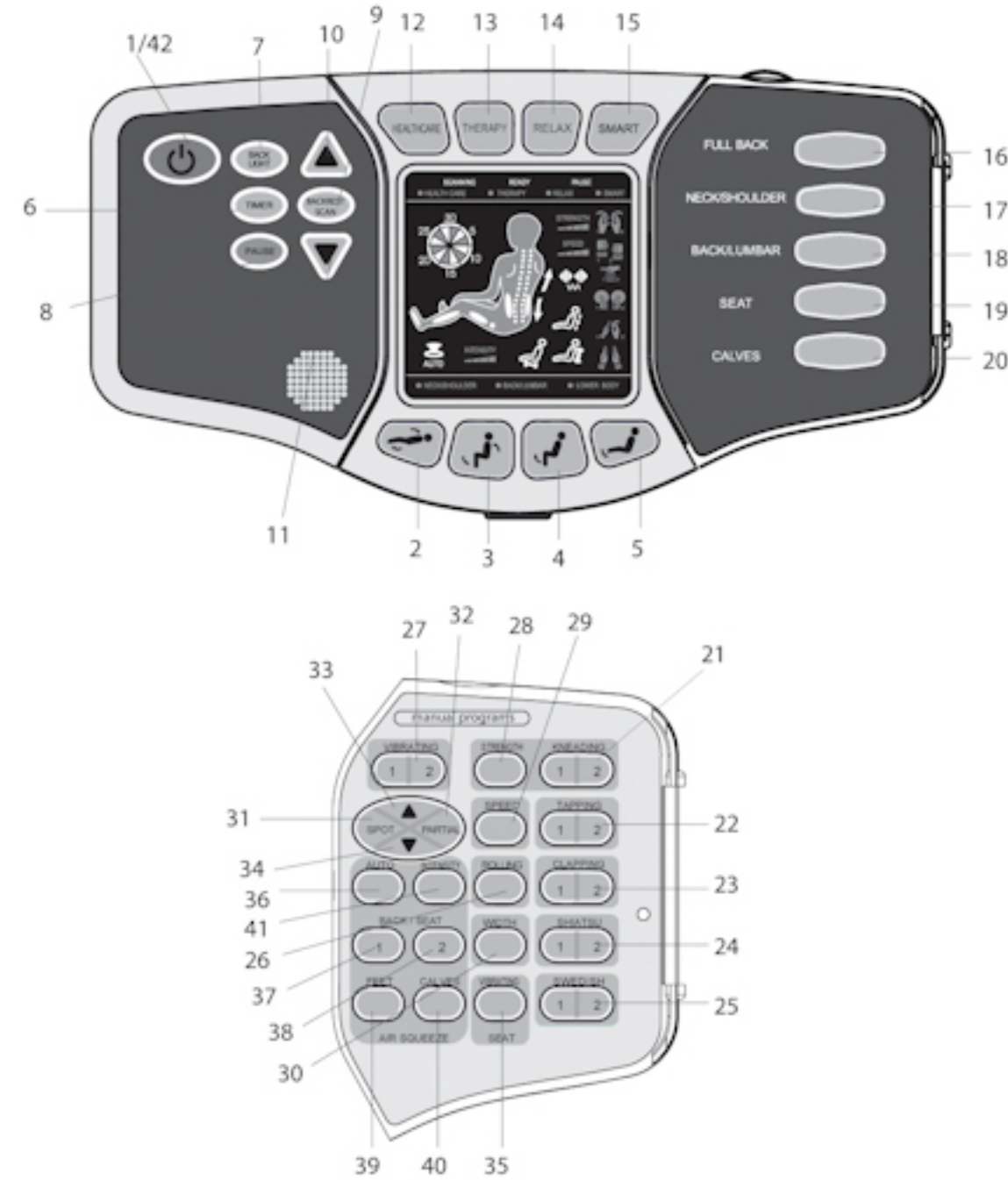Insert the power plug into an electrical outlet. Turn the power switch to the "I" position for idle mode. The following functions are available in idle mode: ON/OFF, Backlight, Recline, Sit Up, Lift Leg, Recline Leg.

Massage chair control buttons

1. ON/OFF

Start/Stop Control Button

Press this button and the backlight of the remote controller will blink to indicate that the massage chair is in idle mode. The backrest will recline, and the footrest will rise and extend to its maximum length. When the chair has stopped moving, any massage function button may be pressed. Press the ON/OFF button again to stop all massage functions and to return the massage rollers to their original position.

Note: Press any function button to begin the massage functions, when any function button is pressed, the backrest and the footrest adjustments will stop.

2/3. Recline/Sit Up

Backrest and Legrest Angle Adjustment Control Buttons

Press one of these buttons once to adjust the backrest and the legrest to a comfortable position. Press one of these buttons again to confirm the position.

Note: Press any function button to begin the massage functions, and the recline and sit up adjustments will stop.

4/5. Lift Leg/Recline Leg

Legrest Angle Adjustment Control Button

Press this button once to adjust the legrest to a comfortable position. Press the button again to confirm the position. Press any control button to adjust recline, sit up, lift leg or recline leg.

Note: 1) Press any function button to begin massage functions and the legrest stop lifting and reclining meanwhile the legrest lifts and reclines. 2) Avoid damaging the footrest, the legrest retracts to its minimum length and extends when it touches ground.

6. Timer

Time Interval Control Button

When a massage function is activated the time is automatically set for 15 minutes. Every time this button is pressed, the time interval will change in the following sequence: 20>25>30>5>10>15 minutes.

Note: The backlight of the remote controller will automatically light up followed by five beeps when one minute remains on the timer.

7. Backlight

Backlight ON/OFF Control Button

Press this button to turn on and off the controller backlight. The backlight automatically turns off after four minutes.

8. Pause

Pause Control Button

Press this button to stop all massage functions and to pause the chair to idle mode. Press the pause button again to resume the original massage function.

9. Backrest Scan

Auto Scan Control Button

Press this button to activate the automatic scan that will position the massage rollers according to the user's height. When the massage rollers are positioned, a beep alerts the user that the manual backrest Up/Down adjustment buttons may be pressed within five seconds to further adjust the massage rollers.

Note: User should sit with shoulders flat against the back of chair for an accurate scan.

10/11. ▲/▼

Manual Backrest Up/Down Control Buttons

Press these buttons to manually adjust the massage rollers.

Full Body Auto Program

12. Healthcare

13. Therapy

14. Relax

15. Smart

Press this button for a brief demonstration of the features of the CZ-610 Massage Chair.

Note: The auto scan function is automatically activated when an auto program is selected. When an auto program is in progress, the back and seat vibration function can be activated and the air squeeze intensity can be adjusted.

Localized Auto Program

A precise massage on a specific part of the body.

16. Full Back Start/Stop Control Button

17. Neck & Shoulder Start/Stop Control Button

18. Back & Lumbar Start/Stop Control Button

19. Seat Start/Stop Control Button

20. Calves Start/Stop Control Button

Note: FULL BACK, NECK/SHOULDER and BACK/LUMBAR auto programs can be combined separately with SEAT,CALVES auto programs. However, the NECK/SHOULDER auto program can't be combined with FULL BACK and BACK/LUMBAR auto programs.

Manual Back Massage

21. Kneading

Kneading Function Start/Stop Control Button – 2 modes (KNEADING 1/KNEADING 2);

4 levels of speed

Press this button once to activate the kneading function for a full back massage at medium intensity.

Press the button again to stop the kneading massage function. Press the KNEADING 1 or the KNEADING 2 button to select the mode desired.

22. Tapping

Tapping Function Start/Stop Control Button – 2 modes (TAPPING 1/TAPPING 2);

3 levels of width (narrow, medium, wide); 4 levels of speed

Press this button once to activate the tapping function for a full back massage at medium intensity.

Press the button again to stop the tapping massage function.

Press the TAPPING 1 or the TAPPING 2 button to select the mode desired.

23. Clapping

Clapping Function Start/Stop Control Button – 2 modes (CLAPPING 1/CLAPPING 2);

3 levels of width (narrow, medium, wide); 4 levels of speed

Press this button once to activate the clapping function for a full back massage at medium intensity.

Press the button again to stop the clapping massage function.

Press the CLAPPING 1 or the CLAPPING 2 button to select the mode desired.

24. Shiatsu

Shiatsu Function Start/Stop Control Button – 2 modes (SHIATSU 1/SHIATSU 2);

3 levels of width (narrow, medium, wide); 4 levels of speed

Press this button once to activate the Shiatsu function for a full back massage at medium intensity.

Press the button again to stop the Shiatsu massage function.

Press the SHIATSU 1 or the SHIATSU 2 button to select the mode desired.

Note: Kneading, tapping, clapping and Shiatsu may be combined.

25. Swedish (light kneading & light tapping)

Swedish Function Start/Stop Control Button – 2 modes (SWEDISH 1/SWEDISH2)

Press this button once to activate the Swedish function for a full back massage.

Press the button again to stop the Swedish massage function. Press the SWEDISH 1 or the SWEDISH 2 button to select the mode desired.

26. Rolling

Rolling Function Start/Stop Control Button – 3 levels of width (narrow, medium, wide)

Press this button once to activate the rolling function for a full back massage.

Press the button again to stop the rolling massage function.

Note: Tapping, clapping, Shiatsu, Swedish and spot functions cannot be combined with the rolling function.

27. Vibrating

Vibration Function Start/Stop Control Button – 2 modes (VIBRATING 1/VIBRATING 2)

Press this button once to activate the vibration function for a full back massage.

Press the button again to stop the vibration massage function.

Press the VIBRATING 1 or the VIBRATING 2 button to select the mode desired.

28. Strength

Kneading Speed Adjustment Control Button

Press this button to adjust the speed when the kneading function is activated.

29. Speed

Speed Adjustment Control Button

Press this button to adjust the speed when the tapping, clapping and Shiatsu functions are activated.

30. Width

Width Adjustment Control Button

Press this button to adjust the width when the tapping, clapping, Shiatsu and rolling functions are activated.

Spot or Partial Back Massage

31. Spot

Spot Back Massage

Press this button once to activate the massage rollers.

Press the ▲ or ▼ button to select a spot massage for any part of the back.

Press the button again to stop the spot massage function and to resume the full back massage function.

Note: The spot massage is not operable when the rolling function is activated.

32. Partial

Partial Back Massage

Press this button once to activate the massage rollers. The massage rollers perform the partial massage with an 8 cm (3.1") range of movement.

Press the ▲ or ▼ button to select a partial massage for any part of the back.

Press the button again to stop the partial massage function and to resume the full back massage function.

33/34. ▲/▼

Position Adjustment Buttons for Spot or Partial Massage

When the spot or partial massage functions are activated these buttons will move the massage rollers.

Press and hold the ▲ or ▼ button to move the massage rollers up or down.

Release the or button when the desired position is achieved.

Setting Seat Vibration

35. Vibrating

Seat Vibration Function Start/Stop Control Button

Press this button once to activate the seat vibration function. Press again to change to another seat vibration mode.

Press this button a third time to stop the seat vibration function.

Setting Air Squeeze

36. Auto

Auto Air Squeeze Control Button

Press this button once to activate the seat, back, feet, and calves auto air squeeze function.

Press the button again to stop the auto air squeeze function.

Note: Press the INTENSITY button to adjust the intensity while the auto air squeeze function is activated.

37/38. Back /Seat

Back and Seat Air Squeeze Start/Stop Control Button –2 modes (AIR SQUEEZE 1/ AIR SQUEEZE 2); 4 levels of intensity

Press this button once to activate the back and seat air squeeze functions.

Press the button again to stop the back and seat air squeeze functions.

Press the AIR SQUEEZE 1 or the AIR SQUEEZE 2 button to select the mode desired.

39. Feet

Feet Air Squeeze Start/Stop Control Button – 4 levels of intensity

Press this button once to activate the feet air squeeze function.

Press the button again to stop the feet air squeeze function.

40. Calves

Calves Air Squeeze Start/Stop Control Button – 4 levels of intensity

Press this button once to activate the calves air squeeze function.

Press the button again to stop the calves air squeeze function.

41. Intensity

Air Squeeze Intensity Adjustment Control Button – 4 levels of intensity

Press this button to adjust the intensity of the air squeeze function while the air squeeze functions are activated.

Note: The vibration and air squeeze functions can be combined with kneading, tapping, clapping, and Shiatsu functions.

42. Stop

Start/Stop Control Button

Press this button to stop all massage functions. The massage rollers will resume their original position.

Switch the power button under the chair to the "O" position and disconnect the power cord plug from the electrical outlet.

Troubleshooting & product maintenance

This massage chair is controlled by a microcomputer and is susceptible to power surges of +10 to -15%.

The following situations may display in the event of a power surge:

Situation	Solution
The massage chair is connected to the power supply but is not functioning.	<ol style="list-style-type: none">1. Check to see that the plug is well connected and working.2. Make sure that the massage chair switch is turned to the "I" position.
The remote control works, but the massage chair is not functioning.	The massage chair may have overheated due to prolonged use. Turn off the power and allow the chair to cool for 30-50 minutes.
Both the remote control and the massage chair are not functioning.	<ol style="list-style-type: none">1. Make sure the chair is plugged into a working outlet.2. Turn the power OFF and then ON again.

If you continue to have problems, please read the warranty card and contact the service department.

Disigned and developed by Anatomico, Italy
Produced in China

Warnings

- Do not use water and other cleaning detergents on the massage chair.
- Do not stack, stab, slit or harshly pull on the surface of the massage chair.
- Gently press the buttons during operation and avoid hitting the buttons with hard objects.
- Pull out the power supply plug after each use.
- Do not use the massage chair for more than 30 minutes at a time.
- DO NOT position any body part in between the gaps of the leg rest when the automatic leg extension feature is in operation as demonstrated PIC 1, in order to prevent pinch to you.



Note: If you are remiss to position any body part in between the gaps of the leg rest, please press PAUSE button to stop and then get your body part out rapidly as demonstrated PIC 2.

DO NOT attempt to press the ON/OFF button to stop, as when we press this button, the footrest will automatically shrink its original position and may cause much tighter pinch.

- If you have to rise from the chair to stop accident, injury or damage to others by you or you want to get your feet and calves out when the automatic leg extension feature is in operation, please press PAUSE button to stop as demonstrated PIC 2, the footrest will gradually deflate in some seconds and then quickly get your foot and calves out and rise.

Note: This way as above mentioned just using in emergencies.

Please bear in mind that always return the footrest to its original length and restore it to its original position, then rise from the chair in order to prevent accident, injury or damage to yourself and chair.

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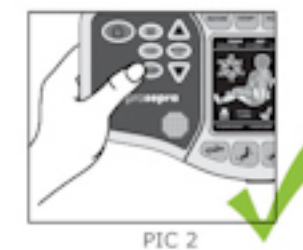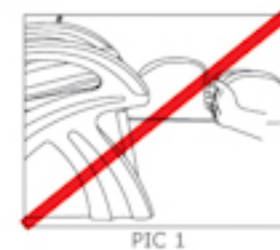Other ANATOMICO products

Massage chair VERDI

Luxurious massage chair VERDI will perfectly add any refined interior, and the fullest range of functions can surprise even the most tempted user. Expensive materials, exclusive design and functionality will ideally underline your individuality and refined taste.

Massage chair PERFETTO

The unique combination of high-tech, design and luxury is inherent in massage chair PERFETTO. The broadest scale of functions, unique technology "Zero gravity", magnificent furnish, a price and quality combination really does PERFETTO by the leader of sales among massage chairs ANATOMICO all over the world.

Massage chair LEONARDO

Massage chair LEONARDO combines high-tech, professional massage and fine design that will allow to look in a new fashion at an interior in your house or office. The convenient and multipurpose mechanism of massage of feet will provide full comfort and pleasure.

Massage chair MARCO

Massage chair MARCO combines high-tech and great functionality. Thanks to popular massage programs and new technologies MARCO is one of leaders in the class. Besides massage of a neck, a back and hips the massage chair allows to do massage of calves. Presence of an air massage provides even more effective relaxation and pleasant rest.